User manual

MANUEL D'UTILISATION
GEBRUIKSAANWIJZING
GEBRAUCHSANWEISUNG
MANUALE DI ISTRUZIONI
MANUAL DE INSTRUCCIONES
INSTRUKCJA OBSŁUGI
NÁVOD K OBSLUZE

L70

Preface	4
1. Your product	5
1.1. Compatibility	5
2. Before use	6
2.1. Intended use and indications	6
2.2. Safety instructions and applicable risks	7
2.3. Symbols on the product	8
3. Using your Neoflex system	9
3.1. Assembly	9
4. Maintenance	10
4.1. Time of maintenance	10
4.2. Maintenance instructions	11
4.3. Expected service life	11
4.4. Reuse	11
4.5. End of use	11
4.6. Warranty	12
5. Technical specifications	12

Preface

Congratulations! You are now owner of a Vermeiren Neoflex cushion!

Thanks for your trust in the products of Vermeiren. To support you on the use and operating options of your product, this manual is offered. If you still have questions after reading this manual, do not hesitate to contact your specialist dealer. They will be glad to help you on this subject.

Important note

To ensure your safety, and to prolong the lifetime of your product, please take good care of it according to the instructions in this manual.

This manual reflects the latest product developments. Vermeiren has the right to implement changes to this type of product without any obligation to adapt or replace similar products previously delivered.

Pictures of the product are used to clarify the instructions. Details of the depicted product may deviate from your product.

Information available

On our website <http://www.vermeiren.com/> you will always find the most recent version of the following information. Please consult this website regularly for possible updates.

Visually impaired people can download the electronic version of this manual and have it read out by means of a text-to-speech software application.

	User manual For user and specialist dealer
	User manual of the wheelchair For user and specialist dealer
	EC declaration of conformity

1. Your product

EN

1. Headrest
2. Backrest
3. Arm pads (optional)
4. Seat
5. Leg rest
6. Footplate

1.1. Compatibility

This product is compatible with the following wheelchairs:

- Inovys II (L70 frame type)

2. Before use

2.1. Intended use and indications

In this paragraph a brief description of the intended use of your product is given. Additionally, relevant warnings are added to the instructions in the other paragraphs. In this way we would like to make you aware of the possible misuse that may appear.

- This product is an accessory to a medical device.
- Indications and contraindications: The Neoflex L70 granule seat system is intended for passive wheelchair users, such as elderly people or people with walking difficulties or no walking abilities (e.g. paralysis, loss of limbs, limb defects or deformations, stiff or damaged joints, heart insufficiencies and poor blood circulation, balance disturbances, cachexia (decrease in muscle), ...). The user has the wheelchair pushed by an attendant.
- The L70 granule seat system is designed to be used in combination with the Inovys II wheelchair with T30 wheel system and frame alterations for the granule seat system.
- Please read all technical details and limits of your Neoflex system in chapter 5..
- The warranty on this product is based on normal use and maintenance as described in this manual. Damage to your product caused by improper use or lack of maintenance will cause the warranty to lapse.

2.2. Safety instructions and applicable risks

CAUTION

Risk of injuries and/or damage

- Please read and follow the instructions in this manual. Otherwise you may get injured, or your product may get damaged.

Keep the following general warnings in mind during use:

- Only use the attachment systems provided by Vermeiren. Do not structurally modify the device to which Neoflex is attached.
- Before each use, check the product and attachments to the wheelchair/ chair to make sure they are undamaged and well secured.
- After each use, check the user's skin in the places of contact. Increasing redness and irritation can evolve into more serious injuries such as pressure ulcers.
- In case of damage, do not attempt to repair the product yourself. Contact your specialist dealer.
- Do not modify your product in any way.
- It is not recommended to use the device when smoking.
- If the packaging of your product has been damaged, (unintentionally) opened, or affected by environmental conditions (moisture, heat,...) upon delivery, check the device integrity of your product. When in doubt, contact your specialist dealer.

Any serious incident [MDR (EU) 2017/745 §2 (65)] that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

For future communication, make sure to keep the LOT number from the packaging.

2.3. Symbols on the product

	Declaration of conformity
	Medical device
	Batch number
	Type designation
	Catalogue number
	Do not iron
	Machine-washable at 40 degrees
	Do not wring
	Mild dry cleaning allowed
	Do not bleach
	Do not tumble dry
	Manufacturer
	Attention: important information
	It is advised to read the manual

3. Using your Neoflex system

⚠ WARNING

Risk of injuries

- First read previous chapters, as well as the user manual of your wheelchair, and inform yourself about the intended use. Do NOT use your Neoflex system unless you have read and fully understood all instructions.
- In case of doubts or questions, do not hesitate to contact your local specialist dealer, your care provider or technical adviser to help you with this.

The design of the Neoflex seat cushions ensures better pressure distribution to avoid pressure-related conditions.



Each section of the L70 system can easily be attached or removed by means of zippers. Depending on the user's needs, the pressure on body contours can be altered per section by adding or reducing the granule filling.

The optional, washable teddy cover can be used for increased hygiene, while other options are available to increase comfort and a sense of safety for the user. Most options are also easily attached/removed without the need for tools.

3.1. Assembly

The L70 system is delivered fully assembled on the wheelchair except for the headrest, the side support pillows of the backrest and seat, and the optional arm pads or leg separator.

- The headrest can be mounted and adjusted according to the instructions in the wheelchair's user manual.
- The pillows are easy to attach by means of the zippers. The label on each pillow indicates the specific place and direction of the pillow (see examples below).

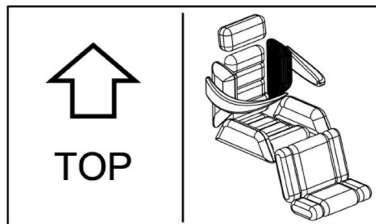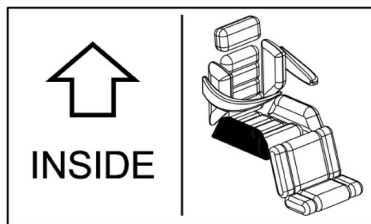4. Maintenance

4.1. Time of maintenance

CAUTION

Risk of injuries and damage

- Repairs and replacements may only be undertaken by trained persons and only genuine replacement parts of Vermeiren should be used.

4.1.1. Before each use

Inspect the following points:

- All parts: present, undamaged and showing no excessive wear, no discolouration.
- All parts: Clean, see § 4.2.1..
- All zippers: closed.

4.1.2. If stored

Make sure that your Neoflex system is stored in a dry environment to prevent mould from growing, or the system from being damaged.

4.2. Maintenance instructions

4.2.1. Cleaning

CAUTION

Risk of damage

- Do not bleach, tumble dry, or iron.

Minor dirt can be cleaned with a cloth. Thorough cleaning by hand wash or a delicate washing programme is advisable depending on the use frequency and the user's needs.

The upholstery can be washed with lukewarm water (max. 40°C) and a mild soap. Do not use abrasive cleaning agents to clean.

The granule seat cushions can be cleaned with a damp cloth (not drenched), and mild household soap if necessary. The teddy covers can be machine washed at 60° C. Do not use abrasive cleaning agents.

4.3. Expected service life

The average service life of your Neoflex system is 5 years. Depending on the frequency of use and maintenance, the service life of your Neoflex product will increase or decrease.

4.4. Reuse

Neoflex products are meant for single patient multiple use only.

4.5. End of use

At end of life, you need to dispose of your Neoflex system according to the local environmental legislation.

4.6. Warranty

The warranty on this product is subject to the general terms and conditions of each country. Please consult you country's website.

5. Technical specifications

The technical details below are only valid for the Neoflex systems at standard settings and optimal ambient conditions. Take these details into account during use. The values are no longer applicable if your system has been modified, damaged, or is severely worn.

Brand	Vermeiren
Product group	Neoflex
Type	L70
Description	Dimensions L70 in combination with Inovys II
Maximum user weight	135 kg
Overall length	1280 mm
Overall width	680 mm 730 mm
Overall height (without headrest)	990 - 1060 mm
Total mass	38 kg
Effective seat depth	460 - 520 mm
Effective seat width	450 mm 500 mm
Backrest height	550 - 600 mm 700 - 750 mm
Footplate to seat distance	360 - 385 - 410 - 435 mm
Use temperature	+5°C - +41°C
Storage and use humidity	30% - 70%
We reserve the right to introduce technical changes. Measurement tolerance +- 15 mm / 1,5 kg / 1,5°	

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**Service registration form**

This product (name):

was inspected (I), serviced (S), repaired (R) or disinfected (D):

By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:
By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:
By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:
By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:
By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:	By (stamp): Kind of work: I / S / R / D Date:

VERMEIREN GROUP

Vermeirenplein 1/15

2920 Kalmthout

Belgium

www.vermeiren.com

Version: G, 2026-09

Basic UDI-DI: 5415174180939Neoflex2H

All rights reserved, including translation.

VERMEIREN
we care **for you**